



All in One

Ready to use self contained swab and diluent system

Swabbing tool "Swabstar" will contribute to the confirmation to effect of the chapter 11, "Cleaning and sanitizing" in FSSC22000, combined "ISO/TS22002-1:2009" and PAS223.



Hygiene & Environmental Testing Swab kit



Unique Swab Head

- ① **Polyester**; Collection efficiency of bacteria compared to the rayon will be up by about 1.5 times. (ELMEX ratio)
- ② **Thin flat design** is effective to sampling for narrow gap places as well as a flat surface area.

With connecting an optional

Extension tool

(patent pend); is effective to sampling for hard-to-reach places, such as deep of machinery or high place.

Stick is able to concatenate 1. 2 or 3 pieces.



Hand and fingers




Knives



For **Allergen** inspection.


Ready to use. Start anytime

Swing stirring after sampling.

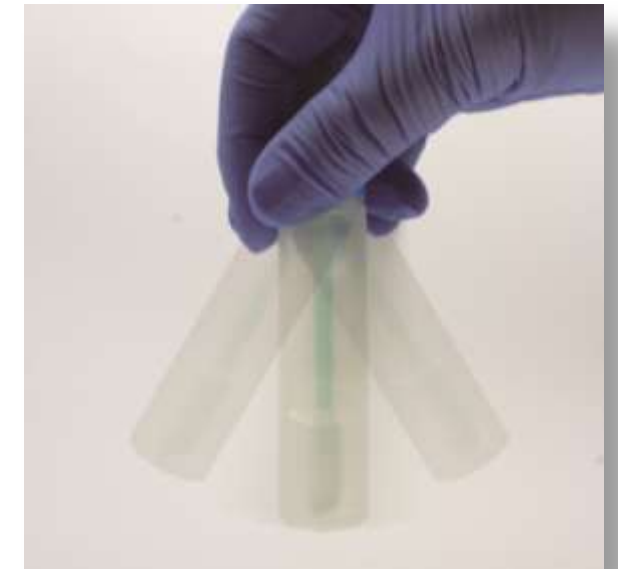

No pipetting require.

Flip open the lid then invert and inoculate media using the measuring gauge on the side by squeezing.

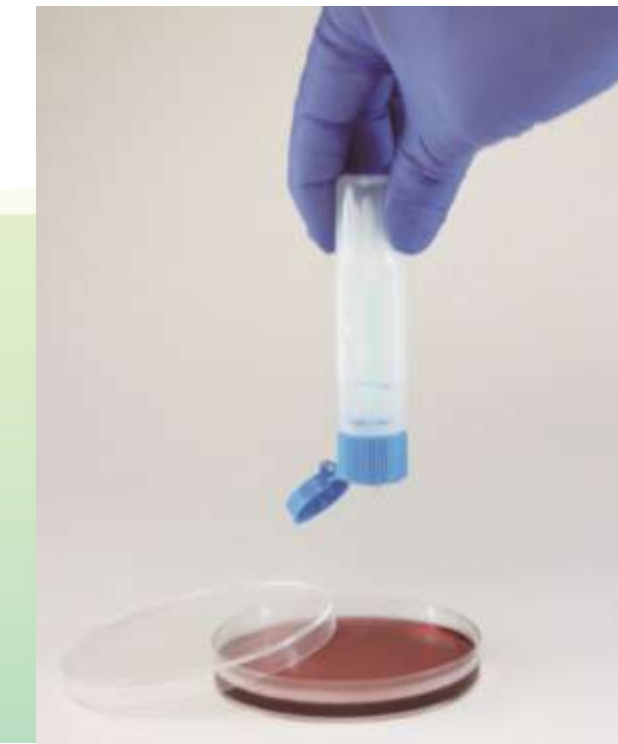



Hygiene & Environmental Testing Swab kit



All in One

Ready to use self contained swab and diluent system



Accessory: Extension tool

Left: capper + stick x 1 Right: capper + stick x 2

Specifications

Code No.	Liquid	Liquid	Pc/box	irradiation	Material				保存条件
					Swab head	Stick	Cap	Bottle	
SS-01-400	Phosphate Buffered Saline (PBS)	10mL	400	Gamma ray	Polyester	Polypropylene	Polypropylene	Polyethylene	2~25°C

- Specifications are subject to be changed for enhancement, without notice.
- This product is for environmental inspection (Swab). Please do not use in any applications.
- When using culture medium put in a Swab-Star container, please sterilize it before discard.
- When discard, please follow your municipality's instructions

Accessory

Code No.	Product name	Pc/pack	Length (cm)	Material
SE-4	Extension tool	4 pairs (Capper x 4pcs, Stick x 4pcs)	Capper + Stick x 1pcs = 155 Capper + Stick x 2pcs = 280 Capper + Stick x 3pcs = 405	ABS

- Specifications are subject to be changed for enhancement, without notice.

Swabbing tool "Swabstar" will contribute to the confirmation to effect of the chapter 11, "Cleaning and sanitizing" in FSSC22000, combined "ISO/TS22002-1:2009" and PAS223.



2-2, Ichigayasadohara-cho, Shinjuku-ku, Tokyo 106-0842
Tel: +81-(0)3-5229-7911 www.elmex.co.jp